




Quadrant Uni-SE-Bond

etching, priming, bonding and desensitising in 1 step

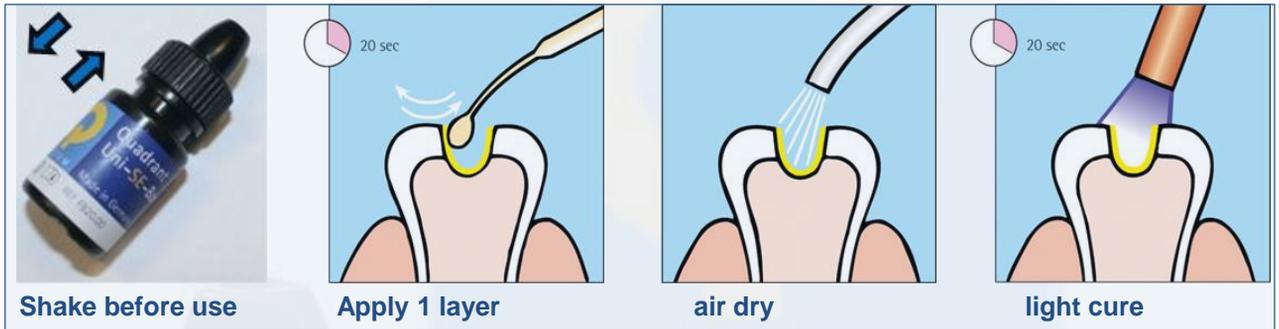

Advantages

- Strong and durable bonding through 4-META and mild etch
- Simply apply 1 layer + air dry + light cure = very simple procedure
- Store at room temperature = easy and reliable

Indications

- Bonding of direct composite restorations
- Bonding of indirect restorations in combination with light-cured adhesive cements
- Sealing of hypersensitive elements

Quadrant Uni-SE-Bond is a 7th generation, light-cured, self-etching bonding system used in combination with light-cured restoration materials. **A unique “built-in” desensitiser** prevents post-operative sensitivity.

Quadrant Uni-SE-Bond permeates deep into the surface, making the hybrid layer thicker and resulting in **improved bonding**. The homogeneity of the hybrid layer is further strengthened (even on vertical walls) by a nanofiller, resulting in less sensitivity **and edge crack-free restorations**.

